

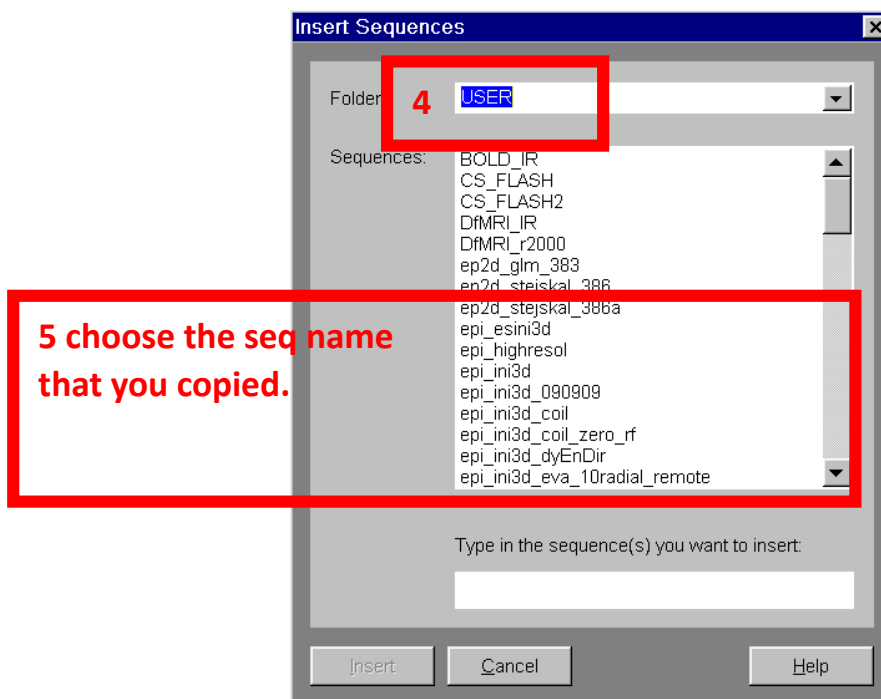


Sequence import (binary files are in dir: "seq")

- 1) Copy 2 binary files under the dir: *C:/MedCom/MriCustomer/seq/*
- 2) insert the sequence

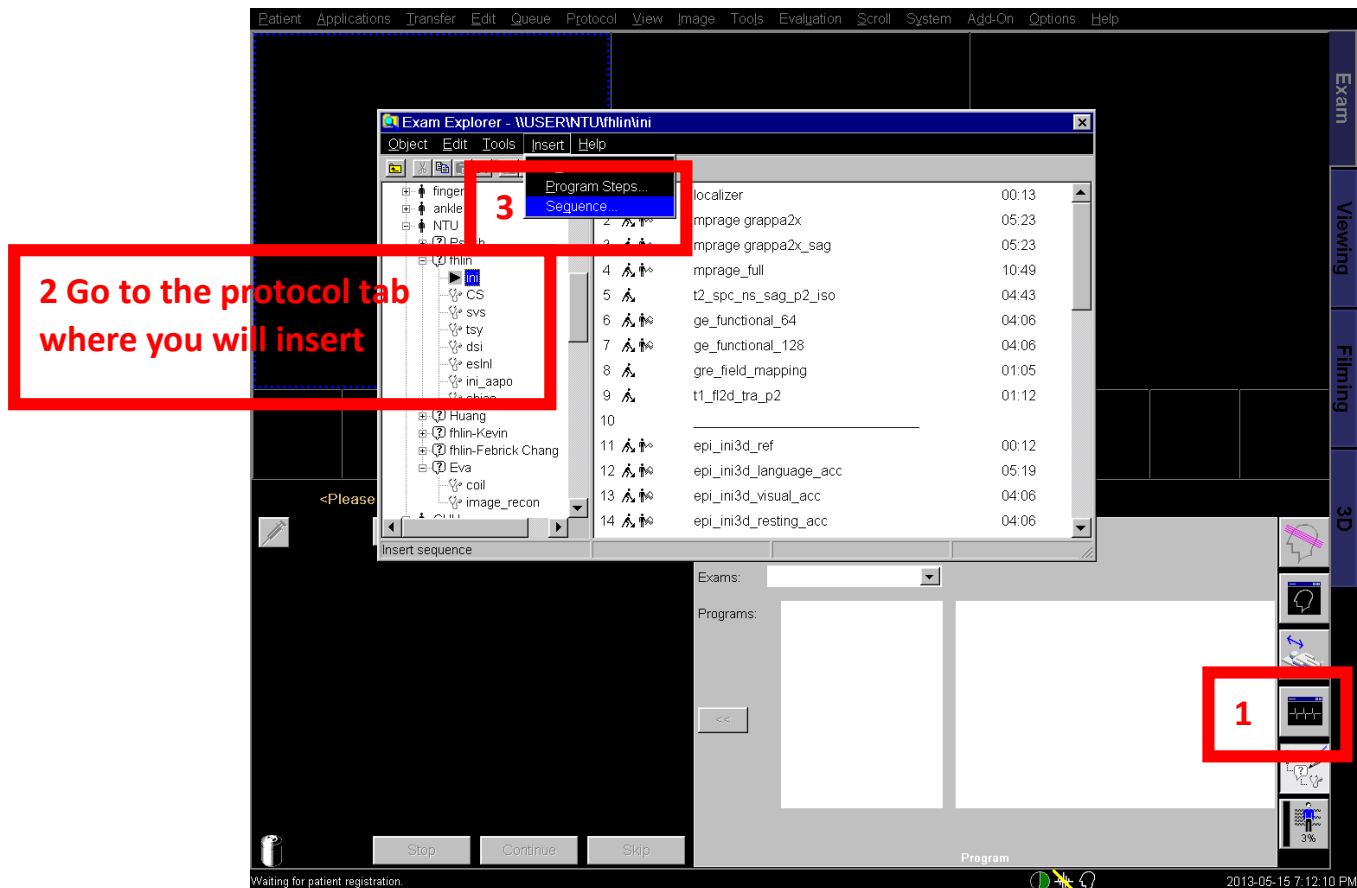

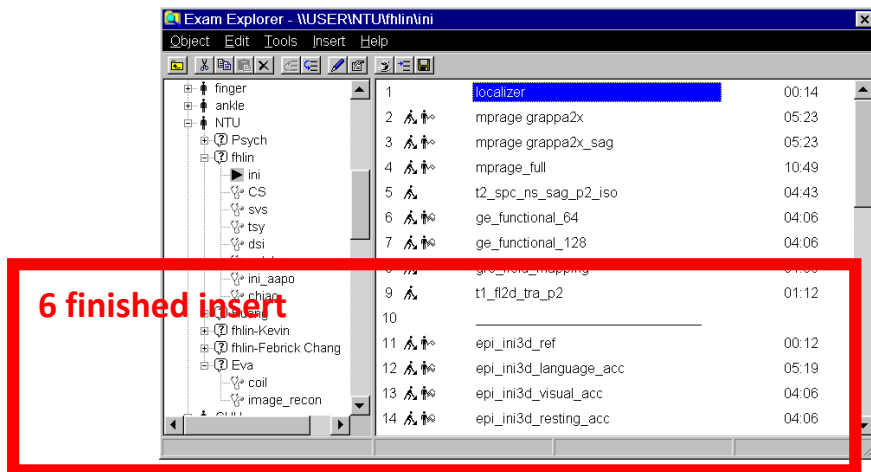

Sequence protocols:

Copy the *seq (smsini)* in protocol tab and rename (one for ref and one for acc)

1) Reference scan (2 segments, takes few sec, size ~ 200MB)

1. Slice = 20
2. Phase enc. Dir = R>>L
3. Fov read = 210 mm
4. TR = 50 ms
TE = 27.5 ms
5. Flip angle = 25 deg
6. Measurements = 10
7. Base resolution = 42
8. Bandwidth = 3970 Hz/px
9. Reference Scan box checked

2) Acc scan (2 segments, takes few mins, size ~ XX GB)

1. Slice = 20
2. Phase enc. Dir = R>>L
3. Fov read = 210 mm
4. TR = 50 ms (the whole brain sampling rate is 10Hz = 0.1 sec),
TE = 27.5 ms
5. Flip angle = 25 deg
6. Measurements = 2460 (= 60 dummy scan + 2400 scan = 246 sec for example)
7. Base resolution = 42
8. Bandwidth = 3970 Hz/px
9. Reference Scan box unchecked

To achieve higher scan number, please use "Averages" option

10. Averages = 1,2,3,.....

For example: total 20 min = $20 \times 60 / (0.1 \text{ sec}) = 12000$ scans

set "Measurements = 2000" and "Averages = 6",

let "Measurements" * "Averages" = total scan

Reconstruction (matlab code files are in dir: “*recon*”):

Using matlab toolbox that we provide

```
>> [ ref, acc ] = MBrecon ( 'meas_ref.dat','meas_acc.dat', 0.0001 );
```